




## + Sterile Preparation Mats

Formatted for Aseptic Pharmacy  
AMD Sterile Triple Wrapped Consumables

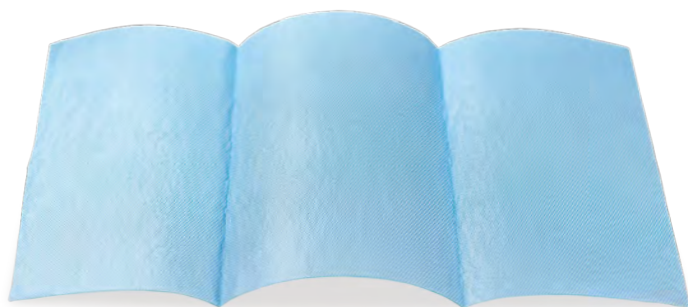

260mm x 200mm



### + ORDERING INFORMATION

Sterile Preparation Mat	
Description	Pack of 1
Sterile Preparation Spill Mat (260mm x 200mm)	Z01A01A
Sterile Preparation Spill Mat, full sized (270mm x 410mm)	Z01A02A

There are two mat sizes available:

- Small - 260mm x 200mm - folded into one third
- Large - 270mm x 410mm - folded in half

### + PRODUCT INFORMATION

The AMD cleanroom sterile preparation mats have following properties:

- Suitable for use in up to ISO 5 areas as it is very low linting

- Latex free
- High liquid holding capacity, porous top layer, bonded reservoir layer and impermeable back layer

The mats are folded into our pouches, one mat per pouch, and each pouch is double outer wrapped.



## + Sterile Kwill (Luer Slip)

Formatted for Aseptic Pharmacy  
AMD Sterile Triple Wrapped Consumables

### + ORDERING INFORMATION

Sterile Kwill (Luer Slip)		Quantity per Pouch	
Description	Length	5	11
Kwill (Luer Slip)	127mm	Z05A01K	Z11A01K



### + PRODUCT INFORMATION

Sterile Luer Slip Kwill for drawing up of liquids, syrups and other viscous solutions, when attached to a syringe. It is virtually transparent and will not kink when in use.

The Kwill is 127mm in length, with an external diameter of 4.95mm and bore diameter of 3.94mm.

Each Kwill is individually packed within the primary pouch, this is to reduce the risk of contamination by a gloved hand whilst removing one from an open pouch.

### Aseptic Medical Devices

Newmarket Drive, Derby DE24 8SW | +44 (0) 1332 755622  
sales@asepticmedical.com | www.asepticmedical.com  
Aseptic Medical Devices is a trading name of Riverside Medical Packaging Company Limited

